


Specialty Growth Hormone Medications at CHOC

Genotropin ™ MiniQuick


Click for Additional Patient Information



1. Wash hands with soap and water.

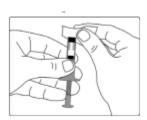


- 2. Open the Genotropin ™ MiniQuick and wipe the rubber stopper with an alcohol pad.
- 3. Peel off the seal from the back of the injection needle, leave the inner and outer needle cover on the needle.
- 4. Attach the needle onto the Genotropin ™ MiniQuick by pushing the back end of the needle down onto the rubber stopper at a 90 degree

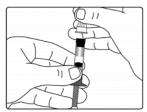

angle. Then, screw the needle clockwise to attach it.

- 5. Hold the Genotropin ™ MiniQuick with the needle pointing straight up and turn the plunger rod clockwise until it does not turn further. This mixes the hormone powder and the liquid.
- 6. Gently tip the MiniQuick in a see-saw motion to dissolve the powder completely (do not shake the pen). Check to see if the liquid is clear and colorless. If it is not clear or colorless, **do not inject the**MiniQuick and call Pfizer Bridge Program (1-800-645-1280).
- 7. Clean the skin at and around the injection site.

- a. Stomach (avoid 2 inches around the navel)
- b. Front of thighs
- 8. With your free hand, gently grasp or pinch a large area of the skin at the injection site and hold firmly for the entire injection procedure.
- 9. Remove the outer and inner needle cap. Position the needle end of

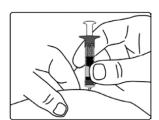

the MiniQuick at a right angle (90 degrees) to the skin.

- 10. Push the plunger rod in as far as it will go to inject all the medicine. Leave the MiniQuick in the injection site for a few seconds after pushing the plunger down in order to ensure all the medication has been delivered.
- 11. After it has been removed, place the outer needle cap back on the needle, and dispose of the used MiniQuick in a sharps container.

Additional Resources:

Hormone Health Network 1-800-HORMONE (1-800-467-6663) hormone.org

Human Growth Foundation 1-800-451-6434 hgfound.org

The Magic Foundation 1-800-362-4423 magicfoundation.org